

# Product Specification

# SCHOTT

## Nexterion® Glass B ultrasonically cleaned (EU)

Dok-Nr.:	PS-HBM-M-0004
Version:	3
Page:	1/2
Date:	23.04.2009

Product Type:	Nexterion® Glass B
Cleaning:	Ultrasonically cleaned
Item Number:	1055132

### 1.1 Description


Nexterion® Glass B is a "microscope slide" sized piece of laser cut borosilicate glass. This glass is manufactured from high purity ingredients using the microfloat process to produce extremely flat sheets. The borosilicate glass exhibits a high chemical resistance, as well as extremely low background fluorescence in the Cy3™ and Cy5™ emission wavelengths of 570 nm and 670 nm respectively. The slides are ultrasonically cleaned under alkaline conditions, after which each slide is checked for conformity. The approved slides are packed into plastic boxes.

### 1.2 Physical Specifications

Nexterion® Glass B is provided with laser cut edges.

Item	Specification
Thickness [mm]	1.00 mm ± 0.050 mm
Length [mm]	75.60 mm ± 0.10 mm
Width [mm]	25.0 mm ± 0.10 mm
Flatness* [mm]	≤ 0.050 mm
Roughness	typical values below 1 nm (RMS)
Transmittance	in the visible light range: > 90 %

\*Flatness = summary effect of thickness deviation (intra-slide) and warp (concave, convex).

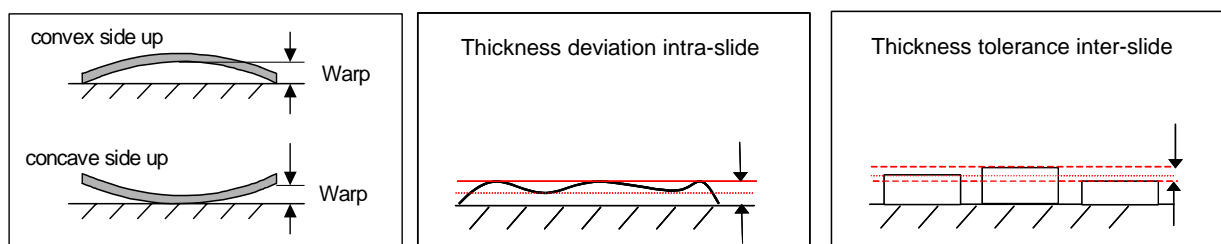

	generated by	reviewed by	reviewed by	reviewed by	authorized by
Name	Ulrike Faessler	Susann Schneider	Stephanie Walther		Christian Jabschinsky
Date		23.04.2009	23.04.2009		23.04.2009

# Product Specification

# SCHOTT

## Nexterion® Glass B ultrasonically cleaned (EU)

Dok-Nr.:	PS-HBM-M-0004
Version:	3
Page:	2/2
Date:	23.04.2009

### 1.3 Surface Specification

Item	Specification	Admissible count
Cleanliness	< 25 surface defects ( $\geq 20 \mu\text{m}$ ) visual inspection under oblique illumination	
Elongated defects (e.g. scratches)	Width $\leq 50 \mu\text{m}$ ; Length $\leq 2 \text{ mm}$	Ignore
	Width $> 50 \mu\text{m}$ ; Length $> 2 \text{ mm}$	None
Spot type defects (e.g. open bubbles, digs, chips, dross, drip, tos, bos)	$< 100 \mu\text{m}$	Ignore
	$> 100 \mu\text{m}$	None
Edge chips	Length $\leq 100 \mu\text{m}$ ; Width $\leq 100 \mu\text{m}$	Ignore
	Length $> 100 \mu\text{m}$ ; Width $> 100 \mu\text{m}$	None
Crack	Visible	None

### 1.4 Contamination

Contaminations like haze, free cullets, dust, foreign particles are permissible, if they are removable by using a standard washing processes for thin glass (alkaline detergent) at the customer's site. Cleanable defects will be ignored and are permissible without any size restriction.

#### For technical assistance please contact:

Europe / Asia – Pacific  
SCHOTT Technical Glass Solutions GmbH  
Otto-Schott-Straße 13  
07745 Jena  
Germany  
Phone: +49-(0)3641-681-4069  
Fax: +49-(0)3641-681-4970  
E-Mail: [coatedsubstrate@schott.com](mailto:coatedsubstrate@schott.com)  
[www.schott.com/nexterion](http://www.schott.com/nexterion)